

# Cooler Insert

## Refrigerate Vaccine at 2°C to 8°C



### Receiving & Unpacking Steps

- 1. Immediately** open cooler and follow TagAlert instructions below.
- 2. Compare** packing slip to cooler contents to confirm they match. (See sample below.)
- 3. Inspect** vaccine and report any tampering (e.g., broken seal). See back to unpack.
- 4. Place vaccine** in appropriate storage immediately upon receipt.
- 5. Repack** cooler for refurbishing following instructions on back.
- 6. Return** cooler via UPS.



### TagAlert Instructions

- Remove TagAlert temperature monitor from box. TagAlert is inside inner box next to vaccine.
- Follow instructions on blue card:
  - Stop TagAlert Indicator.** Press and hold **Start & Stop** button until Stop icon appears.
  - Read Indicator Status.** Determine if appropriate temperature was maintained during transit.
  - TagAlert Monitor Concerns?** Please call **1-877-836-7123**.



INSTRUCTIONS FOR CONFIRMING PROPER SHIPPING TEMPERATURES

#### Receiver Instructions

##### 1. Stop TagAlert® Indicator

Press and hold the START & STOP button for 5 seconds until Stop icon appears.

##### 2. Read Indicator Status

✓ **Accept:** Left arrow points to Green check mark. Ready for Use.

✗ **Reject:** Right arrow points to Red X and numbers 1, 2 and/or 3 are visible—a temperature alarm has triggered. **QUARANTINE PRODUCT** and call +1-877-836-7123 or your state/local immunization program **right away**.

T170852 Rev. A



### Review Packing Slip

NDC		Customer P.O.	Material Description Manufacturer	MFR Lot#	Exp. Date	VFC Doses	317 Doses	State Doses	CHIP Doses	Order Qty	Ship Qty	Unit Price	Extended Price
19515-0810-52 5017406		394261	FLU; SYR; 10-pack (TIV) PED CDC GLAXO SMITHKLINE	AE2J7	06/16/2025	10	0	0	0	10	10	\$15.22	\$152.20
<b>Total</b>											10	10	\$152.20

Verify **MFR Lot#** and **Ship Qty** on packing slip against cooler contents.  
**Ship Qty** is the total number of vaccine cartons.

# Cooler Insert Refrigerate Vaccine at 2°C to 8°C



Unpack & Return All EcoFlex Coolers

Do Not  
Dispose



**Step 1:** Open cooler and remove top gel packs.




**Step 2:** Remove products from inner box.



**Step 3:** Ensure inner box is back in cooler.



**Step 4:** Place gel packs back on top.

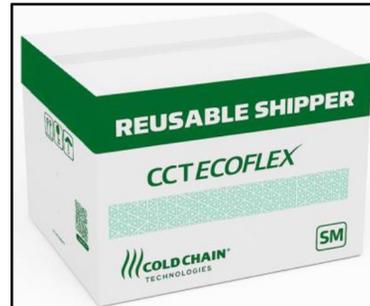

**Step 5:** Fold outer flaps. Fold in remaining flaps with return shipping label exposed. Tape seam.



**Step 6:** UPS automatically picks up within 1-2 days. Return all EcoFlex coolers.



1-800-377-4877



**Do Not Return  
KoolTemp Coolers**

**Please Dispose**